

# PLEGO



**VERMEIREN**  
we care for you

# Easy to fold



# 24,5 kg

Total weight



# Explore the world with Plego



Powered wheelchairs can sometimes be hard to store due to their bulky size or their large total mass, especially when you're travelling. Introducing Plego, a foldable lightweight powered wheelchair that can be folded easily. Even when folded, the wheelchair stays upright and can be manoeuvred not unlike a suitcase.

Whether you are just out and about for a short trip to the shops or exploring a city, Plego can be easily used to explore whilst extending your mobility. The 2 x 250W brushless geared motors are not only quiet and powerful when driving, they are very low maintenance and durable in use. In need of a little extra urban comfort?

Plego has a variety of options such as extra storage or lights.



## MADE IN EUROPE

**READY TO GO  
WHENEVER YOU WANT!**

Built upon 65+ years of Belgian expertise and entirely crafted in Europe, this powered wheelchair embodies our commitment to the utmost standards of quality, safety, and comfort.

Discover more about Plego's features and benefits to boost your life!



Easy  
to Store

In just a few steps, Plego transforms into a compact and easy to store wheelchair. The complete tool free system can be folded by a single person without any strain or without the hassle of removing the battery.

Resembling a travel suitcase, Plego has four additional wheels keeping the chair upright even when folded. They also allow you to maneuver the chair to its designated storage space.

Although Plego's motors pack a punch, the total weight of the wheelchair has been reduced to a mere 24.5 Kg. Road trips have never been easier. Fold into a compact size, lift into the boot of a car without any strain and you are ready to go! Our designers made sure that this wheelchair can be taken wherever it is needed.



# Up to 21 km\* 6 km/h

With a maximum speed of 6 km/h, Plego is a perfect choice to go on adventures with your friends and family. It has a comfortable driving speed to keep up whilst having a conversation comfortably during a walk. The lightweight Lithium-ion battery gives you a driving range of up to 21 km\*. Perfect for a day out!

## PLEGO POWER

The motor selection for Plego emphasizes quality, durability, and efficiency. The brushless gear motor with 2 x 250W power includes a built-in gear mechanism that offers excellent efficiency and longevity compared to brush geared motors. The low-maintenance motors are not only cost-effective due to their design, but are also discreet and quiet whilst driving.




\*17,5 AH BATTERY

# Extra battery needed?

Ready for a day out? Check your battery status easily whilst seated thanks to our battery status indicator located underneath your seat. The Lithium-ion battery can be locked safely with an anti-theft key or charged separately from your chair. No need to bring your entire chair closer to the outlet. Just take your battery and charge it wherever is most convenient. It is important to note that Plego is approved for air transport according to IATA standards\*

To ease your mind about the driving range or your battery charge, extend your battery life with an additional battery (12,8Ah or 17,5Ah) that can easily be stored underneath your seat with the additional storage bag. Check out the options on page 28 to see how you can improve your Plego.

\*Please consult our manual for the most up-to-date information. Always contact your airline as each company often has different rules.



Removable  
lithium-ion  
battery



### LESS IS MORE

Plego's controller is simple, yet efficient. We chose a quality powered wheelchair controller made by Penny & Giles to offer the best steering comfort possible. The standard settings of the controller ensure a smooth and comfortable driving experience. We made sure that anybody could use Plego. The further the joystick is pushed, the faster the power wheelchair will go. Releasing the joystick automatically applies the brakes and stops the wheelchair. With just a few buttons and easy-to-understand LED lights, Plego can be driven and controlled easily. Just get in and go!

### USER FRIENDLY DESIGN

The controller can be mounted interchangeably left or right and has an optional swing-away feature. Additional optional LED lights for higher visibility can be added and are independent from the controller.



## Adventure time?

No need to stress about turning short corners or navigating narrow shopping aisles. Plego's compactness and small turning radius allow you to drive around without worry when going about your day. The lightweight frame and easy to use folding system are perfect to bring Plego wherever you like to explore urban and indoor areas without worry. No car? We've made sure that Plego can easily stand in a bus or any other form of transportation on its own. Just sit back and relax.



GO FURTHER

We care for  
you

**VERMEIREN**, manufacturer of wheelchairs since 1957, is a major player in the world of medical equipment. Over the years, the VERMEIREN group has specialized in producing wheelchairs, nursing beds, shell chairs, comfort chairs and any other product that may accompany a patient or user in maintaining his autonomy and mobility.

GO ELECTRIC

# At home

When and where you use Plego is completely up to you. This lightweight and foldable wheelchair can be stored when you don't need it around or can be used for all day use in indoor spaces. Navigate narrow spaces at home, on public transport or in urban places with ease.

With a width of only 59,5cm Plego fits most European doorways easily and avoids scratching doorposts or other narrow spaces. No need for your Plego at home? Store your powered chair away dust free with our optional Plego cover.



Plego has been designed with a legacy of more than 65+ years of experience in wheelchairs. Even our standard features provide optimal comfort and versatility for all users. Both the backrest and seat are upholstered in 3D Net fabric to optimize airflow, comfort and breathability.



## THE PERFECT COMPANION FOR ADVENTURE

On the back of Plego you'll find a standard included pocket, perfect for a laptop or some reading material.

Not enough storage? We offer two additional pockets located underneath the seat for more storage options. You'll find that Plego has a strong aluminium frame to withstand daily wear and tear like no other.

The foldable monobloc footplate can be adjusted in two height positions and serves a dual function. Not only does it provide support when seated, it also protects frame and battery when transporting the wheelchair in a folded position.



## C77 Glossy Black

Black, the epitome of timeless sophistication, evokes a sense of mystery and elegance.

It embodies power and authority, often symbolizing formality and prestige. In fashion and design, it acts as a versatile canvas, enhancing the impact of other colors while adding a touch of sleekness to any context.



## C88G Glossy Gold

The radiant allure of gold is synonymous with wealth, luxury, and achievement.

It embodies prosperity and success, illuminating the path to opulence. As a symbol of affluence, gold transcends time and culture, leaving a shimmering legacy of extravagance and grandeur.



## C56G Metallic Blue

Metallic blue represents a fusion of modernity and innovation, embodying the sleek, high-tech aesthetics of the future.

This cool and futuristic color captivates with its sense of precision and efficiency. In various industries, it reflects cutting-edge design and technology, symbolizing progress and adaptability.



## C52 Glossy Turquoise

The calming shade of turquoise embraces tranquility and balance.

Its soothing influence invokes the serene qualities of nature, resembling the pristine waters of tropical paradises. As a color of renewal and rejuvenation, it fosters feelings of well-being and spiritual growth.



# C87G

## Glossy Melon Orange

Melon orange radiates warmth and enthusiasm, infusing surroundings with a burst of energy and joy.

This vibrant hue is associated with creativity and playfulness, making it a perfect choice for lively gatherings and artistic endeavors. It embodies the spirit of celebration and cheerfulness.



# C4G

## Glossy Red

Red ignites the passions within us, representing love, energy, and excitement.

This bold and vibrant hue captures attention effortlessly, whether in the realm of romance, art, or design. It symbolizes determination and courage, resonating with the fiery essence of the human spirit.



**B64 STORAGE BAG**  
UNDERNEATH THE SEAT



**B63 STORAGE BAG**  
BEHIND THE FOOT REST



**SE19 LIGHTS**



**HEADREST**




**ESE45**  
STORAGE COVER

## Options

To cater to a wide range of preferences and needs, we offered multiple options for you Plego. Store your Plego or add some additional features to suit you and your needs.



**EXTRA BATTERY**  
12,8 AH OR 17,5 AH

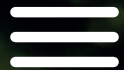

Vermeiren.fm

PLEGO



VERMEIREN  
travel playlist



HEAR ME THROUGH SPOTIFY



# PLEGO

## Technical data

■ Total width	595 mm
■ Total length	1000 mm
■ Seat height	530 mm
■ Seat width	450 mm
■ Seat depth	430 mm
■ Backrest height	460 mm

■ Batteries	12,8 Ah - 24 V or 17,5 Ah - 24 V
■ Motors	2 x 250 W
■ Speed	6 km/h
■ Driving range	15 km or 21 km
■ Weight	24,5 kg
■ Maximum user weight	130 kg

\* Theoretical range calculated in accordance with EU standards under ideal conditions.

\* Speeds and range may vary (approx. 15%).

## VERMEIREN Group

Vermeirensplein 1-15  
B-2920 Kalmthout  
Belgium

[www.vermeiren.com](http://www.vermeiren.com)  
[info@vermeiren.com](mailto:info@vermeiren.com)

For more information, contact your local dealer

